



WellsForm™ Tricuspid Annuloplasty Band



Genesee
BIOMEDICAL, INC.

Shaped by Surgeons.
Built for Better Outcomes.

WellsForm™ Tricuspid Annuloplasty Band

Engineered for Tricuspid Repair — Semi-Rigid and Flexible by Design

Designed For:



- ✓ Tricuspid valve repair procedures
- ✓ Open, port-access, and robotic-assisted approaches
- ✓ Surgeons seeking annular flexibility with segmental control
- ✓ Anatomies requiring precise septal support with anterior adaptability

Clinical Considerations

- ✓ Designed to balance structural support with anatomical adaptability
- ✓ Follows the natural movement of the right atrioventricular junction
- ✓ Ideal for tricuspid repairs involving distorted or dilated annuli
- ✓ Compatible with robotic and minimally invasive surgical workflows

Key Features:

- ✓ Semi-rigid design at the septal and aortic segments supports reshaping the tricuspid annulus
- ✓ Flexible anterior section allows the band to move harmoniously with the free wall of the ventricle
- ✓ Unique curvature conforms to variable planar directions of the tricuspid annulus
- ✓ Braided construction enhances needle penetration and suture security
- ✓ Large radiopaque marker ensures easy fluoroscopic identification
- ✓ Holder includes quick-release cut marks for streamlined band release

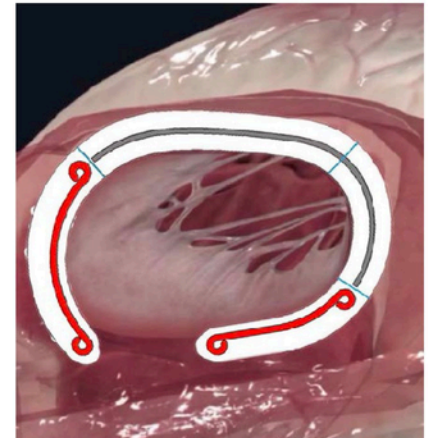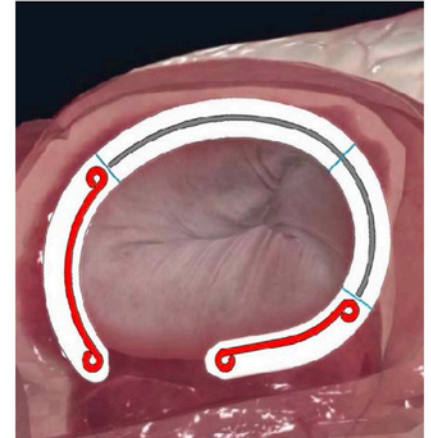


WellsForm™ – Ordering Reference

Ordering Tip: Order using model number WF-##

Example: WF-32 for a WellsForm™ Band, size 32 mm


Catalog #	Product Description
WF-26	WellsForm™ Annuloplasty Band – Size 26 mm
WF-28	WellsForm™ Annuloplasty Band – Size 28 mm
WF-30	WellsForm™ Annuloplasty Band – Size 30 mm
WF-32	WellsForm™ Annuloplasty Band – Size 32 mm
WF-34	WellsForm™ Annuloplasty Band – Size 34 mm
WF-36	WellsForm™ Annuloplasty Band – Size 36 mm

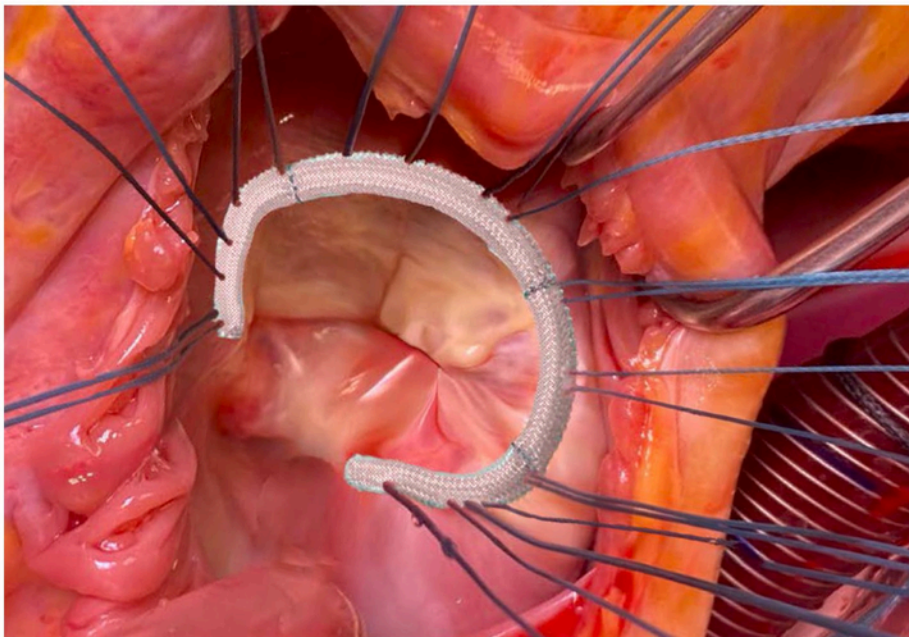

WellsForm™ Accessories – SKU Reference

Note: Accessories are sold separately. Order using the model number.

Catalog #	Accessory Description
WFBS	Set of Reusable Silicone Sizers
GTSS	Set of Stainless-Steel Sizers (includes 2 handles)
GRH-U	Reusable Handle
GRH-UXL	Extra-Long Reusable Handle (35 mm)

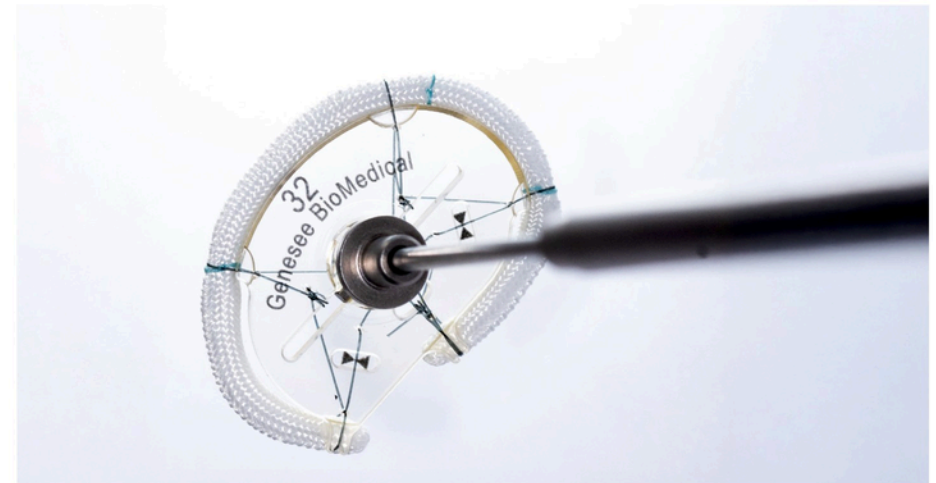

Indications: The WellsForm™ Annuloplasty Bands are indicated for the surgical reconstruction or remodeling of diseased or damaged tricuspid valves. The band provides support for and restricts expansion of the tricuspid annulus.

Warnings: Only surgeons who have received adequate training to determine whether incompetent, stenotic, or diseased heart valves are capable of being repaired or replaced should use this device. Refer to the Instructions for Use (IFU) for a full list of warnings, precautions, and adverse events.

CAUTION: Federal law (USA) restricts this device to sale by, or on the order of, a physician.

Image Credits: Watercolors courtesy of Mr. Francis C. Wells MS, MA (Cantab), FEACTION, F.R.C.S.(Eng) Cambridge, UK.



WellsForm™ Tricuspid Annuloplasty Band | Semi-Rigid Support for Anatomical


Genesee Biomedical



Genesee
BIOMEDICAL, INC.

WELLSFORM™ TRICUSPID

ANNULOPLASTY BAND

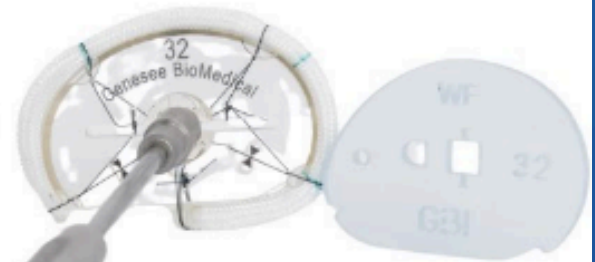

SUBSCRIBE



Watch on YouTube



Annuloplasty Accessories Quick Reference

Simplify. Standardize. Support.

SKU Reference Table

Catalog #	Accessory Description	Compatible Products
FRBS	Set of 9 Flexible Silicone Sizers (Sizes 24–40)	FlexForm™, WellsForm™, TruForm™
GSC	Annuloplasty Scale	TransForm™, FlexForm™, WellsForm™, TruForm™
GSC-2	Set of 2 GSC Annuloplasty Scales and 2 handles (1 GRH and 1 GRH-U)	TransForm™, FlexForm™, WellsForm™, TruForm™
GSS-26/28	Stainless Steel Sizer, double-ended 26/28	FlexForm™, WellsForm™, TruForm™
GSS-30/32	Stainless Steel Sizer, double-ended 30/32	FlexForm™, WellsForm™, TruForm™
GSS-34/36	Stainless Steel Sizer, double-ended 34/36	FlexForm™, WellsForm™, TruForm™
GSS-38/40	Stainless Steel Sizer, double-ended 38/40	FlexForm™, WellsForm™, TruForm™
GSS-4	4 Stainless-Steel Sizers (26/28, 30/32, 34/36, and 38/40) with 2 GRH-U handles	FlexForm™, WellsForm™, TruForm™
GRH	Tricuspid Stainless Steel Sizers, double-ended	FlexForm™, WellsForm™
GRH-U	Reusable Handle	TransForm™, FlexForm™, WellsForm™, TruForm™
GRH-UXL	Extra-Long Handle (35 cm)	TransForm™, FlexForm™, WellsForm™, TruForm™
TFRS	Set of Reusable Silicone Sizers	TransForm™
TFSS	Set of Stainless-Steel Sizers (includes 1 GRH and 1 GRH-U)	TransForm™
WFBS	Flexible Silicone Sizers, set of 6, Sizes 26-36	FlexForm™, WellsForm™

Genesee BioMedical offers reusable sizing tools and handles to support precise valve repair across our annuloplasty portfolio. While most accessories are shared among FlexForm™, WellsForm™, and TruForm™, **TransForm™ utilizes its own sizers** to ensure accurate mitral remodeling.*

Refer to individual product sections or Instructions for Use (IFU) for full cleaning and compatibility details.

All Genesee annuloplasty accessories featured in this section are fully **reusable** and designed for repeat sterilization

* Accessories must be resterilized before each use. Refer to the product's Instructions for Use (IFU) for full warnings, precautions, and reprocessing instructions.
CAUTION: Federal law (USA) restricts these devices to sale by, or on the order of, a physician.

Contact

Shaped by Surgeons.
Built for Better Outcomes.



Genesee BioMedical, Inc.
700 W. Mississippi Avenue, Unit D5
Denver, CO 80223-4509



Phone:
303-777-3000



Website:
www.geneseebiomedical.com



General Inquiries:
info@geneseebiomedical.com

“Innovation begins with listening. Every great product we build starts with a surgeon’s challenge.”
— Matt Dambeck, CEO Genesee BioMedical, Inc

